

# Author Manuscript

## Faculty of Biology and Medicine Publication

**This paper has been peer-reviewed but does not include the final publisher proof-corrections or journal pagination.**

Published in final edited form as:

**Title:** Preparation of a blood culture pellet for rapid bacterial identification and antibiotic susceptibility testing.

**Authors:** Croxatto A, Prod'hom G, Durussel C, Greub G

**Journal:** Journal of visualized experiments : JoVE

**Year:** 2014 Oct 15

**Issue:** 92

**Pages:** e51985

**DOI:** 10.3791/51985

In the absence of a copyright statement, users should assume that standard copyright protection applies, unless the article contains an explicit statement to the contrary. In case of doubt, contact the journal publisher to verify the copyright status of an article.

**Title:**

Preparation of a blood culture pellet for rapid bacterial identification and antibiotic susceptibility testing

**Authors:**

Croxatto Antony, Bugnon 48, Institute of Microbiology, University Hospital Center and University of Lausanne, 1011 Lausanne, Switzerland. [antony.croxatto@chuv.ch](mailto:antony.croxatto@chuv.ch)

Prod'hom Guy, Bugnon 48, Institute of Microbiology, University Hospital Center and University of Lausanne, 1011 Lausanne, Switzerland. [guy.prod'hom@chuv.ch](mailto:guy.prod'hom@chuv.ch)

Durussel Christian, Bugnon 48, Institute of Microbiology, University Hospital Center and University of Lausanne, 1011 Lausanne, Switzerland. [christian.durussel@chuv.ch](mailto:christian.durussel@chuv.ch)

Corresponding author: Greub Gilbert, Bugnon 48, Institute of Microbiology, University Hospital Center and University of Lausanne, 1011 Lausanne, Switzerland.  
[gilbert.greub@chuv.ch](mailto:gilbert.greub@chuv.ch), Tel : +41213144979, fax : +41213144060

**Keywords**

Blood culture, bacteriology, identification, antibiotic susceptibility testing, MALDI-TOF MS.

**Short Abstract**

A rapid bacterial pellet preparation from a positive blood culture can be used as a sample for applications such as identification by MALDI-TOF, Gram staining, antibiotic susceptibility testing and PCR-based test. The results can be rapidly communicated to clinicians to improve the outcome of patients suffering from bloodstream infections.

## Long Abstract

Bloodstream infections and sepsis are a major cause of morbidity and mortality. The successful outcome of patients suffering from bacteremia depends on a rapid identification of the infectious agent to guide optimal antibiotic treatment. The analysis of Gram stains from positive blood culture can be rapidly conducted and already significantly impact the antibiotic regimen. However, the accurate identification of the infectious agent is still required to establish the optimal targeted treatment. We present here a simple and fast bacterial pellet preparation from a positive blood culture that can be used as a sample for several essential downstream applications such as identification by MALDI-TOF MS, antibiotic susceptibility testing (AST) by disc diffusion assay or automated AST systems and by automated PCR-based diagnostic testing. The performance of these different identification and AST systems applied directly on the blood culture bacterial pellets is very similar to the performance normally obtained from isolated colonies grown on agar plates. Compared to conventional approaches, the rapid acquisition of a bacterial pellet significantly reduces the time to report both identification and AST. Thus, following blood culture positivity, identification by MALDI-TOF can be reported within less than an hour whereas results of AST by automated AST systems or disc diffusion assays within 8 to 18 hours, respectively. Similarly, the results of a rapid PCR-based assay can be communicated to the clinicians less than 2 hours following the report of a bacteremia. Together, these results demonstrate that the rapid preparation of a blood culture bacterial pellet has a significant impact on the identification and AST turnaround time and thus on the successful outcome of patients suffering from bloodstream infections.

## Introduction

Bloodstream infections and sepsis in hospitalized patients are a major cause of morbidity and mortality. Thus, mortality related to bloodstream infections is observed in about 14% to 37 % of hospitalized patient and may increase to 35% in intensive care units patients <sup>1-3</sup>. The rapid identification of the infectious agent is pivotal to guide optimal antimicrobial treatment and to increase the successful outcome of antimicrobial therapy. The rapid analysis of Gram stains from positive blood culture has already a significant impact on the adaptation of antimicrobial therapy <sup>4,5</sup> but accurate identification of the infectious agent is required to provide the best adapted antibiotic treatment to the patients. For instance, different antibiotic treatment regimens have to be implemented following bacteremia with enterococci and streptococci that are difficult to distinguish by Gram staining. Similarly, identification at the species level is required to detect Gram negative enterobacteria encoding a chromosomal *ampC* gene which confer an increased resistance to  $\beta$ -lactams <sup>6</sup>.

With a positive blood culture, the conventional diagnostic approach is to subculture the infectious agent on different agar plates, which requires several hours of additional incubation prior identification with various approaches including biochemical tests, growth on different selective media and automated microbial identification systems. The time to results of a conventional diagnostic approach is of about 1 to 3 days.

The emergence of the matrix-assisted laser desorption/ionization time-of-flight mass spectrometry (MALDI-TOF) technology for rapid identification of microorganisms has provided a new tool to quickly identify microorganisms from colonies grown on agar plates but also directly from positive blood cultures (Figure 1) <sup>7-10</sup>. The use of MALDI-TOF to identify an infectious agent from blood cultures has significantly reduced the time to results to a few minutes instead of the hours and days required by traditional methods. As discussed by Croxatto et al. <sup>11</sup>, the efficiency of MALDI-TOF identification relies on different parameters including the microorganism's purity and quantity. These two criteria are easily obtained from discrete colonies grown on agar plates but required a pre-analytical treatment for bacterial enrichment and purification from complex samples such as blood culture, which contain multiple cellular and protein components that may interfere with MALDI-TOF identification.

Various microorganisms' isolation methods from blood culture have been used in a number of studies including saponin or other mild detergents method for bacterial extraction <sup>7,12</sup>, serum separator method <sup>8</sup>, lysis centrifugation methods <sup>10</sup> and commercially solutions such as the Sepsityper kit (Bruker Daltonics, Bremen, Germany). Our bacteriology diagnostic laboratory has developed a simple blood-culture bacterial pellet preparation based on ammonium chloride erythrocyte-lysis which allow fast identification of bacteria and yeast by MALDI-TOF and automated identification systems (Figure 2) <sup>13</sup>. This blood-culture pellet preparation also provide a sample for other direct downstream applications such as Gram staining, automated PCR-based diagnostic tests such as POCT-PCRs for the rapid detection of methicillin-resistant *Staphylococcus aureus* (MRSA), and antibiotic susceptibility testing with automated AST systems and/or by disk diffusion assays on agar plates (Figure 3).

In this work, we describe the different steps for the preparation of the blood-culture bacterial pellet (Figure 4). We will also describe the protocols for three of the main applications that can be performed on the blood culture pellet: Identification by MALDI-TOF <sup>13</sup>, identification (ID) and antibiotic susceptibility testing (AST) with the automated systems <sup>14</sup> for *Enterobacteriaceae* and staphylococci and automated PCR-based diagnostic test for the detection of MRSA <sup>15</sup>.

## **Protocols**

This protocol has been developed and validated following the research and development processes and ethical rules of our institution before being implemented as a routine tool.

### **1. Preparation of a blood culture bacterial pellet by ammonium chloride erythrocyte-lysing procedure**

#### **1.1. Preparation of a positive blood culture for subsequent centrifugation.**

1.1.1. Sterilize the blood culture cap. Add 70 % ethanol on the cap of the bottle and burn it. Do not perform this step under a laminar flow hood.

1.1.2. Move the bottle to a laminar flow hood. Collect 5 ml of the blood culture with a 20 needle gauge syringe

1.1.3. Add the 5 ml of blood culture in a 50 ml Falcon tube containing 45 ml of sterile water (H<sub>2</sub>O) and mix the sample.

#### **1.2. Preparation of a blood culture bacterial pellet by lysis centrifugation.**

1.2.1. Centrifuge the sample at 1000 x g for 10 min at room temperature.

1.2.2. Remove the supernatant (which contains mainly H<sub>2</sub>O and blood cells) with a laboratory vacuum pump.

1.2.3. Resuspend the pellet in 1 ml of home-made ammonium chloride lysing solution (0.15M NH<sub>4</sub>Cl, 1mM KHCO<sub>3</sub>). Break down the pellet by scraping and by vortexing the tube and centrifuge the sample at 140 x g for 10 min at room temperature.

1.2.4. Discard the supernatant which mainly contains red blood cells debris.

1.2.4.1. If the pellet remained hemorrhagic, resuspend the pellet in 2 ml of sterile and distilled H<sub>2</sub>O to further lyse residual red blood cells and to wash the pellet. Centrifuge the sample at 140 x g for 10 min at room temperature and discard the supernatant.

1.2.5. Resuspend the pellet in 200 µl of H<sub>2</sub>O

#### **1.3. Subculture of blood culture on various selective and non selective agar media.**

1.3.1. Collect 1 ml of the blood culture with a 20 needle gauge syringe.

1.3.2. Inoculate 1 drop of blood culture into four quadrants on blood agar, chocolate agar, McConkey agar and Schaedler agar.

1.3.3. Incubate at 37°C the blood agar and chocolate agar plates in a 5% CO<sub>2</sub> atmosphere, the McConkey agar in a normal atmosphere and the Schaedler agar plate in anaerobic conditions.

1.3.4. Follow bacterial growth on the different media.

1.3.5. Check for bacterial purity.

This protocol has been developed and validated following the research and development processes and ethical rules of our institution before being implemented as a routine tool.

## **2. Identification by MALDI-TOF MS**

### **2.1. Direct identification from the blood pellet.**

2.1.1. Transfer 1  $\mu$ l of the blood pellet on a MALDI target plate and let it dry on a 42°C heating platform.

2.1.2. Overlay the sample with 1 $\mu$ l 70% formic acid and let it dry on a 42°C heating platform.

2.1.3. Overlay the sample with 1 $\mu$ l of MALDI matrix (saturated solution of  $\alpha$ -cyano-4-hydroxycinnamic acid in 50% acetonitrile-2.5% trifluoroacetic acid) and let it dry on a 42°C heating platform.

2.1.4. Proceed to MALDI-TOF MS analysis with a MALDI-TOF MS system as described by the manufacturers.

2.1.5. If the identification score is invalid at the species level (for instance with a score < 2), perform a protein extraction of the blood pellet sample.

### **2.2. Identification after protein extraction from the blood culture pellet.**

2.2.1. Mix 20  $\mu$ l of the pellet with 1ml of 70% ethanol and centrifuge at 13'000 x g for 2 min at room temperature.

2.2.2. Discard the supernatant and add 25  $\mu$ l of 70% formic acid and 25  $\mu$ l of 100% acetonitrile to the pellet. Vortex vigorously to resuspend properly the pellet. Difficult to resuspend pellets can be broken down with pipette tips before vortexing.

2.2.3. Centrifuge at 13'000 x g for 2 min at room temperature. Transfer 1  $\mu$ l of the supernatant (which contain extracted proteins) on a MALDI target plate and let it dry on a 42°C heating platform. Proceed as described in 2.1.2

This protocol has been developed and validated following the research and development processes and ethical rules of our institution before being implemented as a routine tool.

## **3. Bacterial identification and antibiotic susceptibility testing with an automated system.**

3.1. Preparation of a bacterial suspension for bacterial identification and AST with an automated system.

3.1.1. Prepare a plastic tube containing 3ml of sterile 0.45 % NaCl solution compatible with the automated system.

3.1.2. Add a sufficient volume of the blood culture pellet in the 0.45 % NaCl solution to obtain a bacterial suspension corresponding to a McFarland turbidity of 0.6 to 0.8 measured with a densitometer machine. The volume of the blood culture bacterial pellet required to

achieve a McFarland turbidity of 0.6 to 0.8 is highly variable but is usually comprised between 50 to 200  $\mu$ l.

3.2. Inoculate the automated Gram-positive (GP) or Gram-negative (GN) cards for identification (this step is not applied if MALDI-TOF MS identification gave a score value >2) and the automated AST cards for antimicrobial susceptibility testing of Gram-positive cocci and Gram-negative bacteria as described by the manufacturer.

3.3. Subculture the bacterial suspension on blood agar (GP and GN) and MacConkey agar (GN) plates to check the purity of the bacterial suspension.

This protocol has been developed and validated following the research and development processes and ethical rules of our institution before being implemented as a routine tool.

#### **4. Antibiotic susceptibility testing by disk diffusion assay as described by the European committee on antimicrobial susceptibility testing (EUCAST, Version 3.0, April 2013, [www.eucast.org](http://www.eucast.org))**

##### **4.1. Preparation of a bacterial suspension to perform antimicrobial susceptibility testing by disk diffusion assay.**

4.1.1. Prepare a glass tube containing 3ml of sterile 0.9 % NaCl solution.

4.1.2. Add a sufficient volume of the blood culture pellet in the 0.9 % NaCl solution to obtain a bacterial suspension corresponding to a McFarland turbidity of 0.5 measured with a densitometer machine. The volume of the blood culture bacterial pellet required to achieve a McFarland turbidity of 0.5 is highly variable but is usually comprised between 50 to 200  $\mu$ l.

4.1.3. Vortex the bacterial suspension

##### **4.2. Inoculation of bacterial suspension onto Mueller-Hinton agar (MH) or Mueller-Hinton agar-fastidious organisms agar (MHF) (MH supplemented with 5% defibrinated horse blood and 20 mg/L of $\beta$ -Nicotinamide adenine dinucleotide).**

4.2.1. Dip a sterile cotton swab in the bacterial suspension and gently remove excess fluid by turning the swab on the inside wall of the glass tube.

4.2.2. Spread the bacterial suspension evenly onto the entire surface of the MH/MHF agar plate by swabbing in three directions or by using an automated plate rotator.

##### **4.3. Application of antimicrobial disks on inoculated agar plates.**

4.3.1. Apply closely the disks on the surface of the dried inoculated agar plates manually or by using a susceptibility disk dispenser.

4.3.2. Incubate the plate at the appropriate temperature and atmosphere as described in the antimicrobial susceptibility testing EUCAST disk diffusion method (Version 3.0, April 2013, [www.eucast.org](http://www.eucast.org)).

This protocol has been developed and validated following the research and development processes and ethical rules of our institution before being implemented as a routine tool.



## 5. Automated PCR-based diagnostic test for the detection of MRSA

5.1. Using a micropipette, transfer 50 µl of the blood culture pellet into the sample reagent vial provided with the MRSA automated PCR-based diagnostic test kit and vortex during 20 seconds

5.2. Using a 1ml micropipette, dispense the sample into the specimen port of the cartridge. Insert the cartridge into the automated PCR system and start the assay as described by the manufacturer.

### Representative results

In the study performed by Prod'hom *et al.*<sup>13</sup>, bacterial pellets obtained by ammonium chloride lysis centrifugation of 122 positive blood culture from 78 patients were analyzed by MALDI-TOF MS. Out of 122 positive blood culture, 95 (77.9%) were correctly identified at the species level and one (0.8%) at the genus level. The remaining 26 (21.3%) blood culture pellets gave no reliable identification by MALDI-TOF. Among those, 21 were gram positive bacteria including 13 streptococci and 5 coagulase-negative staphylococci. Importantly, 10 out of the 13 unidentified streptococci were *Streptococcus pneumoniae* which often possess a difficult-to-lyse capsule and are hardly distinguished from species of the *S. mitis* group using MALDI-TOF MS. Five Gram-negative bacteria were not identified by MALDI-TOF, among which four difficult-to-lyse encapsulated bacterial species including two *Klebsiella pneumoniae* and two *Haemophilus influenzae*. Thus, correct identification by MALDI-TOF was obtained in 78.7 % of blood culture pellets, which represents an efficient performance compared to the 84.1% of identification obtained from conventional identification of bacterial colonies grown on agar plates<sup>9</sup>.

The performance of the automated Vitek 2 system for bacterial identification (ID) and for AST was tested on 278 bacterial pellets of positive blood culture for Gram-positive cocci in cluster and Gram-negative bacteria. These results were compared to conventional results obtained with the Vitek 2 cards from colonies grown on agar plates following subculture of positive blood culture on appropriate media<sup>14</sup>. Direct identification using Vitek 2 cards can be used when identification of the bacterial pellet by MALDI-TOF MS has failed. Compared to a final identification by MALDI-TOF MS, a correct identification with Vitek 2 ID cards was observed in 99% of *Enterobacteriaceae*, 74% of staphylococci, 71% of non-fermentative Gram-negative bacteria and 13% of other Gram-positive cocci. A misidentification was observed in 11% of total cases whereas 8% of bacterial pellets gave no identification. A total of 220 bacterial pellets including 87 *Enterobacteriaceae* and 133 staphylococci were analyzed for AST with Vitek 2 AST GN26 and AST 580 cards<sup>14</sup>. The results were analyzed according to definitions given by the US Food and Drug Administration (FDA) for interpretive agreement results, which defined a false susceptible as a very major error (VME) and a false resistant as a major error (ME). AST obtained from direct inoculation of blood culture pellet compared to inoculation from colonies grown on agar plates showed 0.1% VME and 0.3% ME for *Enterobacteriaceae*, and 0.7% VME and 0.1% ME for staphylococci. Thus the performance of

Vitek 2 AST cards inoculated directly with blood culture pellets are in agreement with the performance criteria accepted by the FDA.

The performance of the PCR-based GeneXpert MRSA test targeting the *spa*, *mecA* and SCC genes was directly applied on *Staphylococcus aureus* blood culture bacterial pellets identified by MALDI-TOF MS<sup>15</sup>. Based on 106 cases, the detection of methicillin resistance with the GeneXpert MRSA exhibited a 99% sensitivity and a 100% specificity. Interestingly, the median time to report the GeneXpert results from blood culture positivity was equal to 201 min (range 100-430) whereas the median time from MALDI-TOF identification to report of the GeneXpert results was equal to 97 min (range 25-250).

## Figures

**Figure 1:** Pre- and post-MALDI-TOF era. Compared to conventional diagnostic approaches which include an overnight subculture of the sample, MALDI-TOF MS identification can be performed within minutes following samples delivery to the diagnostic laboratory.

**Figure 2:** Workflow of blood culture processing from sampling to MALDI-TOF MS identification. When blood cultures are positive, a bacterial pellet is prepared by ammonium chloride lysis centrifugation. This bacterial pellet is then used for direct identification by MALDI-TOF MS. MALDI-TOF MS identification can be obtained within 30 to 60 min following blood culture positivity.

**Figure 3:** Many downstream applications including Gram staining, MALDI-TOF MS identification, AST using automated systems and/or disk diffusion assays and rapid PCR-based testing such as point of care testing PCRs (POCT-PCR) for the detection of MRSA can be performed directly on the blood culture bacterial pellet. Direct testing from the bacterial blood culture pellets significantly decrease the turn-around time to report ID and AST results which is pivotal to improve the outcome of patients suffering from bloodstream infections.

**Figure 4:** The ammonium chloride lysis centrifugation protocol allows the preparation of an enriched and purified bacterial pellet from positive blood culture broths. Gram staining of positive blood culture of *Escherichia coli*, coagulase-negative *Staphylococcus capitis* and *Streptococcus dysgalactiae* performed before and after the preparation of the blood culture bacterial pellet. The classical morphologies of staphylococci in clusters and streptococci in chains are more easily observed in the native blood culture broth than in the blood culture bacterial pellet. Scale bar = 25  $\mu$ m

## Discussion

Compared to conventional positive blood culture diagnostic approaches, the rapid acquisition of a bacterial pellet by using the ammonium chloride lysis centrifugation approach reduce the time to report identification by 16 to 24 hours and the time to report AST by 24 to 48 hours (Figure 1 to 3).

Rapid introduction of appropriate antibiotic therapy is pivotal to improve the outcome of patients suffering from bloodstream infections. Thus, early identification of the infectious agents within an hour followed by a rapid AST obtained in about 8 to 16 hours represents a major impact in the clinical management of septic patients. The rapid reporting of Gram staining has already been associated with a great impact on antimicrobial therapy and thus with a decrease in patient morbidity and mortality<sup>4,5</sup>. The Gram staining of the blood culture bacterial pellet may be done to increase sensitivity when the Gram staining of the blood culture broth is negative or weakly positive (Figure 4). However, Gram staining performed on the native blood culture broth is recommended to observe typical morphologies such as staphylococci in clusters or streptococci in chains. In a recent prospective study performed on 202 cases of bloodstream infection, MALDI-TOF MS identification from blood culture pellets showed an impact in antibiotic therapy in 35.1% of cases whereas Gram stain had only an additional impact in 20% of cases<sup>6</sup>. Interestingly, the maximum difference of impact between Gram staining and MALDI-TOF MS identification reporting on clinical management was observed when bacteria associated with increased antibiotic resistance such as AmpC-producing *Enterobacteriaceae* were identified by MALDI-TOF MS. Moreover, the use of a rapid PCR-based assay such as the GeneXpert MRSA following *S. aureus* identification by MALDI-TOF MS from blood culture pellets allowed to provide the most appropriate antibiotic therapy in less than 4 hours to patients suffering from *S. aureus* bacteremia<sup>15</sup>. When excluding patients with penicillin allergy, the use of the GeneXpert MRSA allowed a significant reduction of anti-MRSA antibiotics misuse such as glycopeptides from 26.1% to 8.1%. Together, these data suggest that rapid sequential reporting of Gram staining, MALDI-TOF MS identification and rapid PCR-based assay from blood culture pellets are significantly improving the antibiotic therapy and thus the clinical management as well as the outcome of patients suffering from bloodstream infections.

MALDI-TOF MS allowed a 78.7% correct identification of blood culture pellets, which is a very good performance compared to a rate of 84.1% identification obtained by conventional MALDI-TOF MS performed on pure colonies grown on agar plates. More than 50% of bloodstream infections are caused by bacterial species such as *Enterobacteriaceae*, *Staphylococcus aureus*, *Pseudomonas aeruginosa* and enterococci. These bacteria are efficiently identified by MALDI-TOF MS which likely positively influence the good performance of MALDI-TOF MS identification performed directly on blood culture bacterial pellets. The low performance for the identification of some bacterial species by MALDI-TOF MS performed on blood culture pellets is mainly due to similar factors than those observed for the conventional identification from pure colonies grown on agar plates. For instance, several microbial species such as Streptococci belonging to the *mitis* group and coagulase-negative staphylococci are closely related and identified with low accuracy by MALDI-TOF MS. Moreover, a particular composition of the bacterial cell wall of Gram-positive bacteria or the presence of a dense capsule observed in bacterial species such as *Klebsiella pneumoniae* significantly decrease the lysis efficiency and thus the performance of MALDI-TOF MS identification. Finally, the low performance of MALDI-TOF MS for the identification of anaerobic bacteria is mainly due to an incomplete and poor representation of these bacteria in the MALDI-TOF database.

The small reduced performance of MALDI-TOF MS identification from blood culture pellets compared to conventional identification from pure colonies may be due to residual blood protein components or to an insufficient amount of bacteria obtained after lysis-centrifugation.

Similarly, the differences observed with the Vitek 2 cards inoculated directly with blood culture pellets compared to those inoculated with isolated colonies may be caused by residual blood culture components such as proteins, blood cells and blood medium compounds that may interfere with bacterial growth in the Vitek 2 cards and may have a significant impact on both identification and AST. Moreover, the AST results obtained by Vitek 2 systems need to be checked by disk diffusion assays also performed directly on blood culture bacterial pellets to validate the results of some antibiotics known to present significant discrepancies compared to conventional approaches. However, direct inoculation of the Vitek 2 cards with blood culture pellets presents an excellent ID and AST performance for both *Enterobacteriaceae* and staphylococci.

Thus, the bacterial pellet obtained by ammonium chloride lysis centrifugation of positive blood cultures has been validated to directly perform several downstream applications allowing a rapid reporting of essential results such as ID and AST in a significantly reduced TAT compared to conventional approaches. Moreover, additional applications such as extended-spectrum  $\beta$ -lactamases (ESBLs) testing have been validated on similar lysis centrifugation of positive blood culture methodologies<sup>16</sup> suggesting that they could also be applied on the bacterial pellet obtained with the protocol described in this work.

## **Acknowledgments**

We thank the technicians of the bacteriology laboratory of the University Hospital Center of Lausanne for their help to implement the techniques in the laboratory.

## **Disclosures**

The authors declare that they have no competing financial interests.

## References


- 1 Alberti, C. *et al.* Epidemiology of sepsis and infection in ICU patients from an international multicentre cohort study. *Intensive care medicine* **28**, 108-121, doi:10.1007/s00134-001-1143-z (2002).
- 2 Angus, D. C. *et al.* Epidemiology of severe sepsis in the United States: analysis of incidence, outcome, and associated costs of care. *Critical care medicine* **29**, 1303-1310 (2001).
- 3 Vincent, J. L. *et al.* Sepsis in European intensive care units: results of the SOAP study. *Critical care medicine* **34**, 344-353 (2006).
- 4 Barenfanger, J. *et al.* Decreased mortality associated with prompt Gram staining of blood cultures. *American journal of clinical pathology* **130**, 870-876, doi:10.1309/AJCPVMDQU2ZJDPBL (2008).
- 5 Munson, E. L., Diekema, D. J., Beekmann, S. E., Chapin, K. C. & Doern, G. V. Detection and treatment of bloodstream infection: laboratory reporting and antimicrobial management. *J Clin Microbiol* **41**, 495-497 (2003).
- 6 Clerc, O. *et al.* Impact of matrix-assisted laser desorption ionization time-of-flight mass spectrometry on the clinical management of patients with Gram-negative bacteremia: a prospective observational study. *Clin Infect Dis* **56**, 1101-1107, doi:10.1093/cid/cis1204 (2013).
- 7 Meex, C. *et al.* Direct identification of bacteria from BacT/ALERT anaerobic positive blood cultures by MALDI-TOF MS: MALDI Sepsityper kit versus an in-house saponin method for bacterial extraction. *Journal of medical microbiology* **61**, 1511-1516, doi:10.1099/jmm.0.044750-0 (2012).
- 8 Moussaoui, W. *et al.* Matrix-assisted laser desorption ionization time-of-flight mass spectrometry identifies 90% of bacteria directly from blood culture vials. *Clin Microbiol Infect* **16**, 1631-1638, doi:10.1111/j.1469-0691.2010.03356.x (2010).
- 9 Seng, P. *et al.* Ongoing revolution in bacteriology: routine identification of bacteria by matrix-assisted laser desorption ionization time-of-flight mass spectrometry. *Clin Infect Dis* **49**, 543-551, doi:10.1086/600885 (2009).
- 10 Stevenson, L. G., Drake, S. K. & Murray, P. R. Rapid identification of bacteria in positive blood culture broths by matrix-assisted laser desorption ionization-time of flight mass spectrometry. *J Clin Microbiol* **48**, 444-447, doi:10.1128/JCM.01541-09 (2010).
- 11 Croxatto, A., Prod'hom, G. & Greub, G. Applications of MALDI-TOF mass spectrometry in clinical diagnostic microbiology. *FEMS Microbiol Rev* **36**, 380-407, doi:10.1111/j.1574-6976.2011.00298.x (2012).
- 12 Chen, J. H. *et al.* Direct bacterial identification in positive blood cultures by use of two commercial matrix-assisted laser desorption ionization-time of flight mass spectrometry systems. *J Clin Microbiol* **51**, 1733-1739, doi:10.1128/JCM.03259-12 (2013).
- 13 Prod'hom, G., Bizzini, A., Durussel, C., Bille, J. & Greub, G. Matrix-assisted laser desorption ionization-time of flight mass spectrometry for direct bacterial identification from positive blood culture pellets. *J Clin Microbiol* **48**, 1481-1483, doi:10.1128/JCM.01780-09 (2010).
- 14 Prod'hom, G., Durussel, C. & Greub, G. A simple blood-culture bacterial pellet preparation for faster accurate direct bacterial identification and antibiotic susceptibility testing with the VITEK 2 system. *Journal of medical microbiology* **62**, 773-777, doi:10.1099/jmm.0.049361-0 (2013).
- 15 Clerc, O. *et al.* Matrix-assisted laser desorption ionization time-of-flight mass spectrometry and PCR-based rapid diagnosis of *Staphylococcus aureus* bacteraemia. *Clin Microbiol Infect*, doi:10.1111/1469-0691.12329 (2013).
- 16 Nordmann, P., Dortet, L. & Poirel, L. Rapid detection of extended-spectrum-beta-lactamase-producing Enterobacteriaceae. *J Clin Microbiol* **50**, 3016-3022, doi:10.1128/JCM.00859-12 (2012).

# Pre-MALDI-TOF era

# Post-MALDI-TOF era




Gram + cocci  
arranged in chains

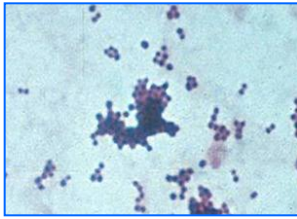


Gram + cocci  
arranged in clusters




Gram - bacilli

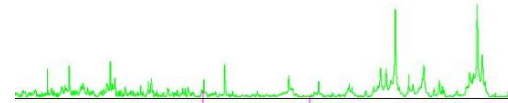

*E. faecium*

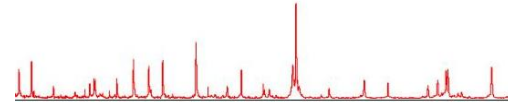

*S. aureus*

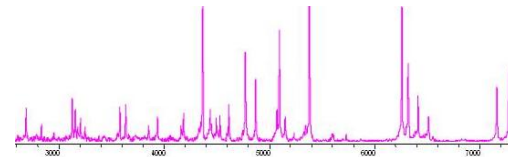


*E. coli*

Overnight  
incubation




*S. aureus*

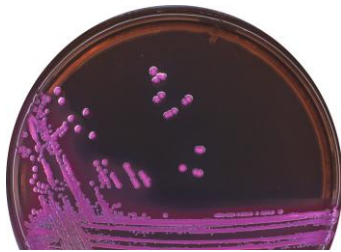

*E. coli*



Day 0  
Day 1

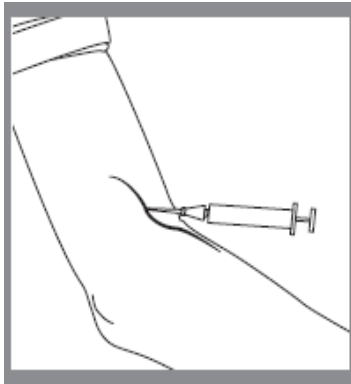

Sampling




Incubation of blood culture bottles



WHEN POSITIVE




Preparation of a bacterial pellet

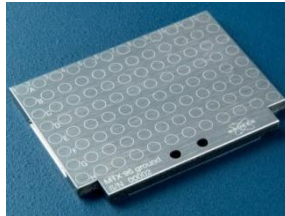

Deposition of bacterial pellet on MALDI microplate



Acquisition of the proteic profile

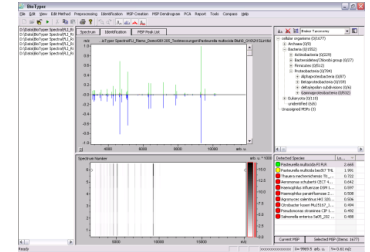


Comparison with a database







Blood culture  
bacterial pellet

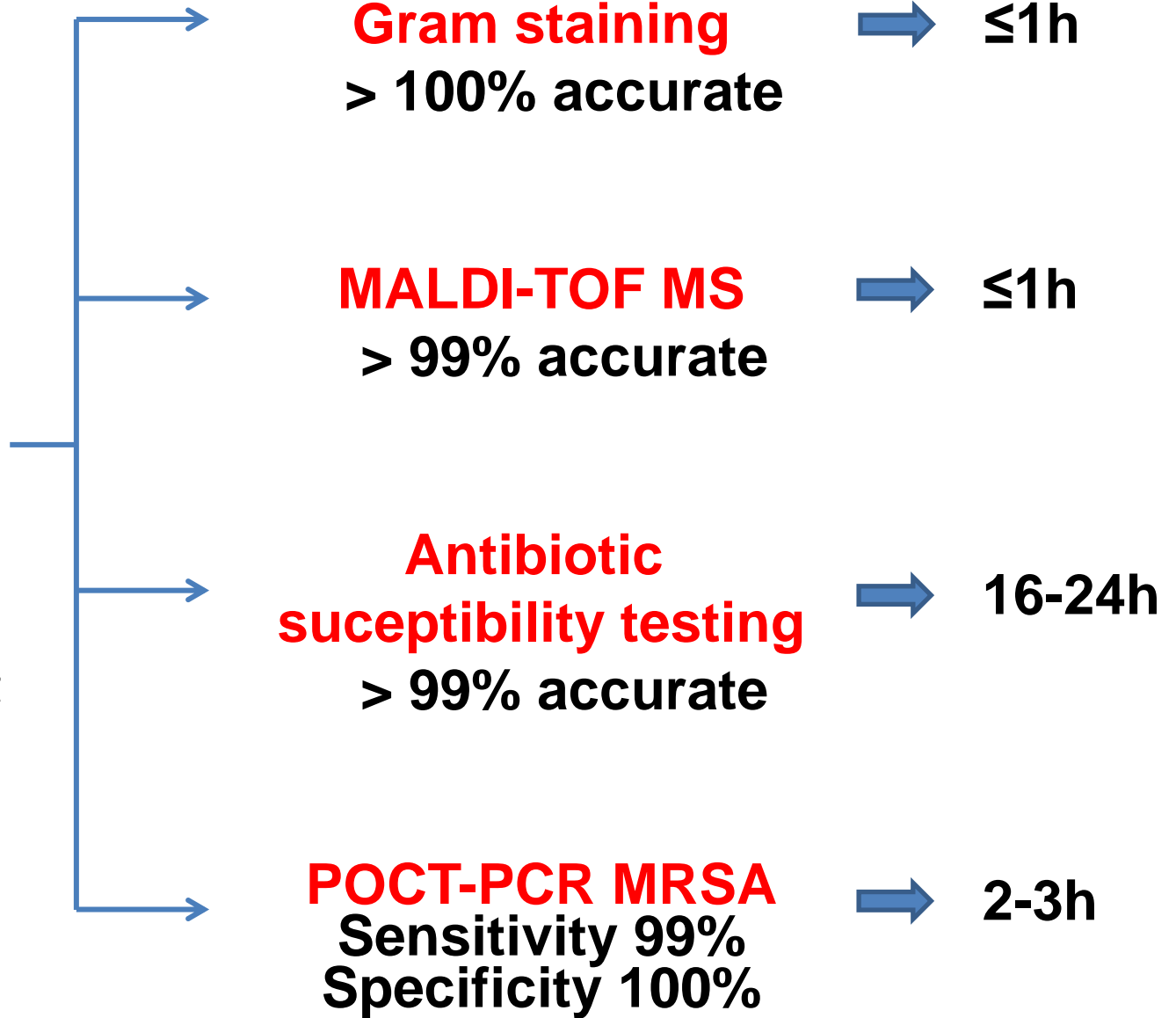



**Blood culture broth**

**Blood culture pellet**



*Escherichia coli*

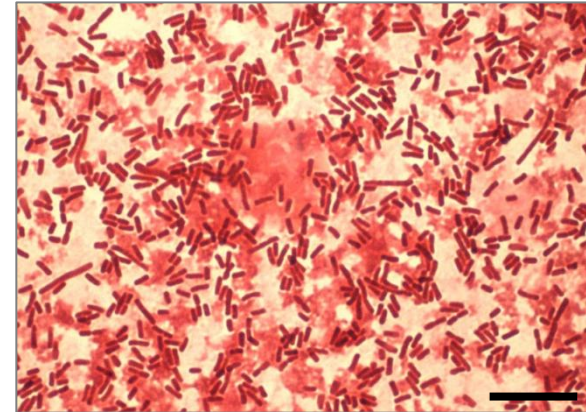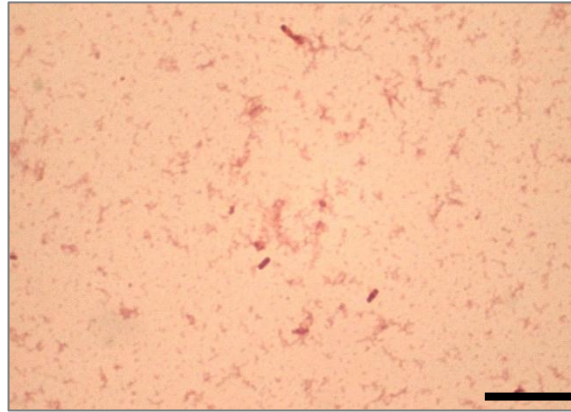

*Staphylococcus capitis*

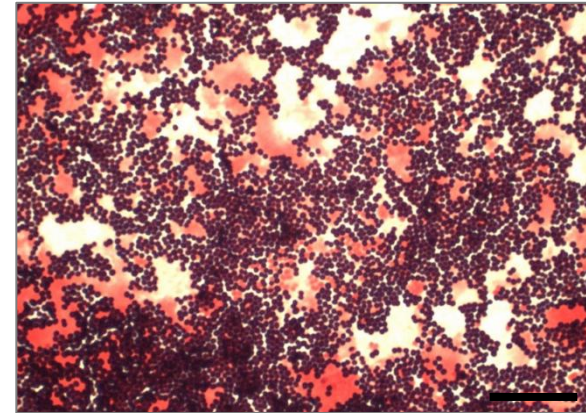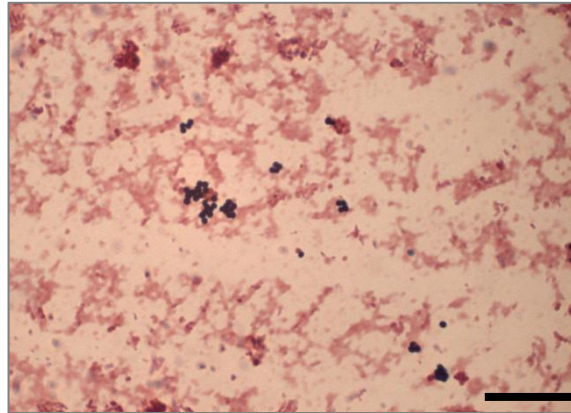

*Streptococcus dysgalactiae*

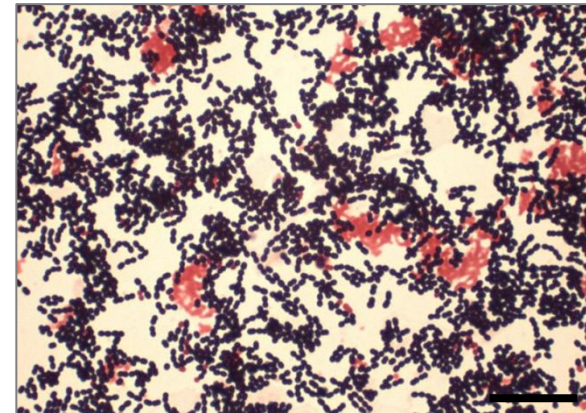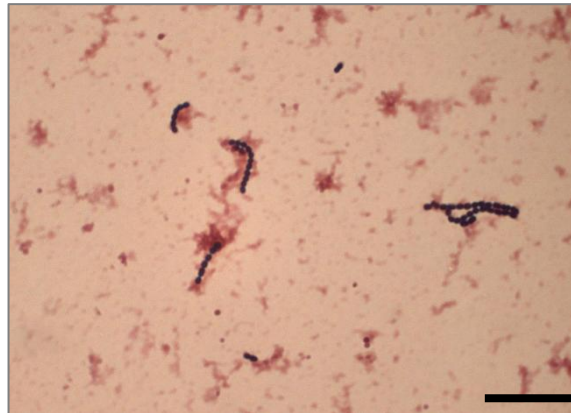

Name of reagent /equipment	Company
20 needle gauge	Terumo, Leuven, Belgium
Ammonium chlorure	Merck, Darmstadt, Germany
Potassium hydrogen carbonate	Fluka, St. Louis, MO
Formic acid	Sigma-Aldrich, St. Louis, MO
$\alpha$ -Cyano-4-hydroxycinnamic acid	Fluka, St. Louis, MO
Acetonitrile	Sigma-Aldrich, St. Louis, MO
Trifluoroacetic acid	Sigma-Aldrich, St. Louis, MO
Vitek 2 60 instrument	Biomérieux, Marcy-l'Etoile, France
Vitek 2 Gram-positive (GP) card	Biomérieux, Marcy-l'Etoile, France
Vitek 2 AST-P580 card	Biomérieux, Marcy-l'Etoile, France
Vitek 2 Gram-negative (GN) card	Biomérieux, Marcy-l'Etoile, France
Vitek 2 AST-N242 card	Biomérieux, Marcy-l'Etoile, France
Xpert MRSA	Cepheid, Sunnyvale, Ca, USA
GeneXpert IV instrument	Cepheid, Sunnyvale, Ca, USA
Microflex LT MALDI-TOF MS instrument	Bruker Daltonics, Bremen, Germany
Densitometer Densicheck instrument	Biomérieux, Marcy-l'Etoile, France

Catalog Number	Comments/Description
NN-2038R	
101145	
60340	
F0507	Flammable, corrosive
70990	Acute toxicity
271004	Flammable, acute toxicity
T6508	Corrosive, acute toxicity
27202	
21342	
22233	
21341	
413391	
GXMRSA-100N-10	
GXIV-4-D	
BDAL microflex LT/SH	
27208	