

# Results SimID.3014001

Priyanshu Bajaj, Markus Kühbach<sup>1</sup>

<sup>1</sup>*Max-Planck-Institut für Eisenforschung GmbH*

## Disclaimer

**This report has been auto-generated using paraprobe-autoreporter.  
It is the responsibility of the author to check the validity and correctness of these results!**

## 1 paraprobe-transcoder

**Here is place for tool-specific comments.**

The dataset contains 184254735 ions in total.


Table 1: Which XML settings were used for paraprobe-transcoder?

Keyword	Value	Unit	Description
Inputfile	R76_30140-v01.pos		
TranscodingMode	1		

## 2 paraprobe-ranger

Table 2: What is the composition of the dataset, i.e. how many ions with particular ion type labelled?

unknown	Sc	Al	Ga	Si	Mg	O:	Xe	AlO:	ScH:	Fe	Co	Ni	AlH:	Total
5839290	519822	177219096	18386	157732	6766	831	5057	6808	80956	9532	1353	2478	386628	184254735

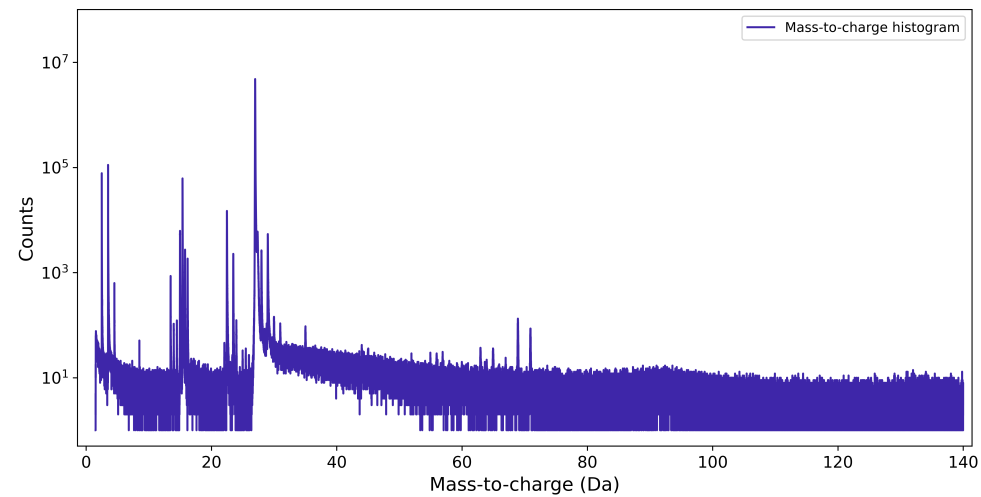

(a)

Figure 1: What is the mass-to-charge histogram/diagram of the entire dataset?

Table 3: Which XML settings were used for paraprobe-ranger?

Keyword	Value	Unit	Description
Inputfile	R76_30140-v01.pos		
TranscodingMode	1		


3 paraprobe-surfacar

4 paraprobe-tessellator

5 paraprobe-spatstat

6 paraprobe-dbscan

PARAPROBE.Surfacr.Results.SimID.3014001.h5  $\alpha$ -shape,  $\alpha = 1.63$

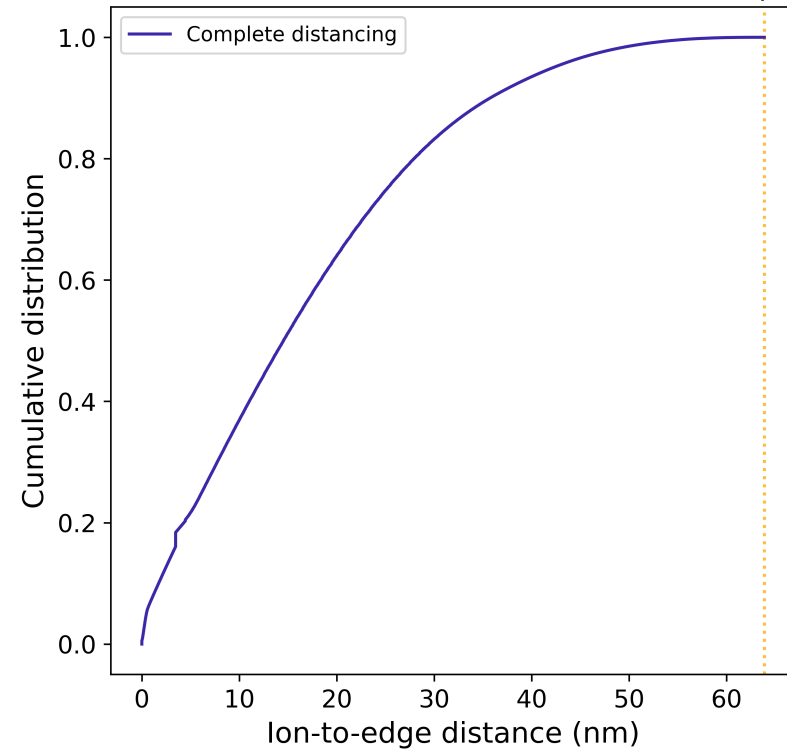

(a)

Figure 2: How many ions are so and so far away from the dataset edge?


## 7 Profiling

Table 4: Which XML settings were used for paraprobe-surfacers?

<b>Keyword</b>	<b>Value</b>	<b>Unit</b>	<b>Description</b>
AdvDistanceBinWidthIncr	1	nm	
AdvDistanceBinWidthMax	1	nm	
AdvDistanceBinWidthMin	1	nm	
AdvIonPruneBinWidthIncr	1	nm	
AdvIonPruneBinWidthMax	1	nm	
AdvIonPruneBinWidthMin	1	nm	
AlphaShapeAlphaValue	0		
DistancingMode	1		
DistancingRadiusMax	0	nm	
InputfileReconstruction	PARAPROBE.Transcoder.Results.SimID.3014001.h5		
RequeryingThreshold	0		
SurfacingMode	1		


Table 5: Which XML settings were used for paraprobe-tessellator?

<b>Keyword</b>	<b>Value</b>	<b>Unit</b>	<b>Description</b>
CellErosionDistance	1	nm	
GuardZoneFactor	5		
IOCellNeighbors	0		
IOCellProfiling	1		
IOCellShape	0		
IOCellVolume	1		
InputfileDistances	PARAPROBE.Surfacers.Results.SimID.3014001.h5		
InputfilePSE	PARAPROBE.PeriodicTableOfElements.xml		
InputfileReconstruction	PARAPROBE.Transcoder.Results.SimID.3014001.h5		
IonsPerBlock	5		
SpatialSplittingMethod	0		

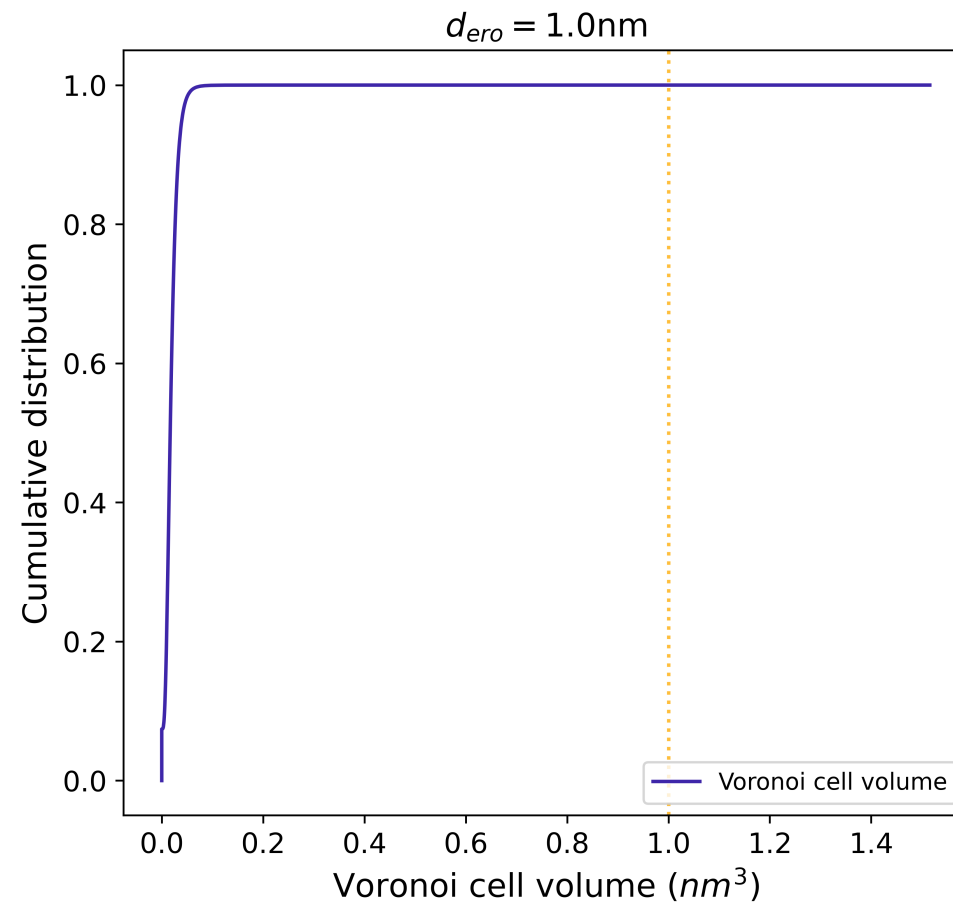


(a)

Figure 3: What is the distribution of volume for the Voronoi cells?

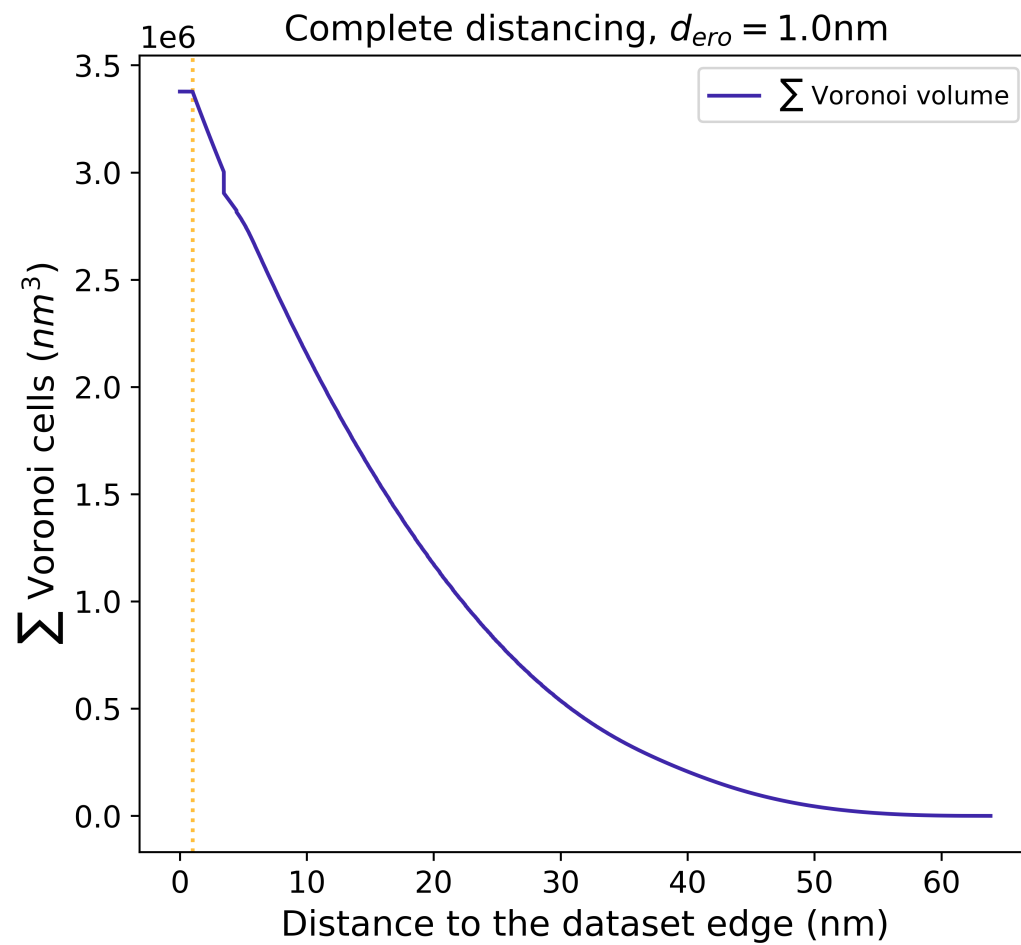


(a)

Figure 4: How much accumulated volume of Voronoi cells remains when successively eroding Voronoi cells from the dataset edge towards the interior?

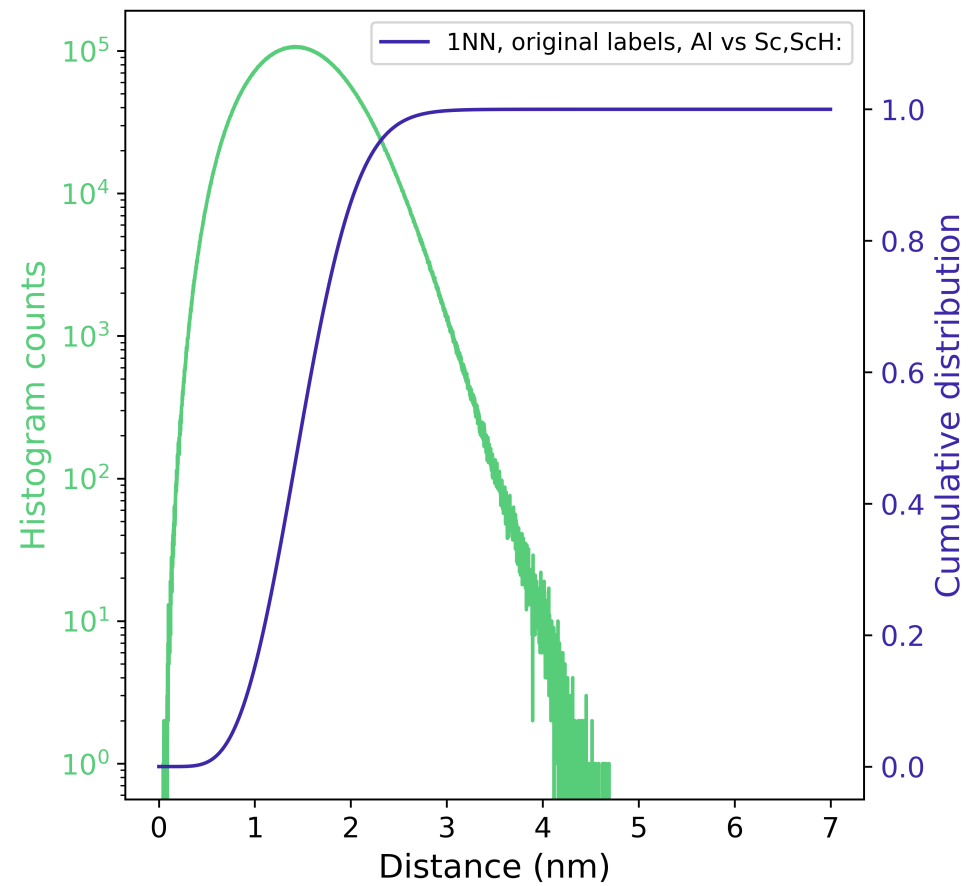

(a)


Figure 5: What is the kth nearest neighbor distribution function for particular ion types about specific ion types?

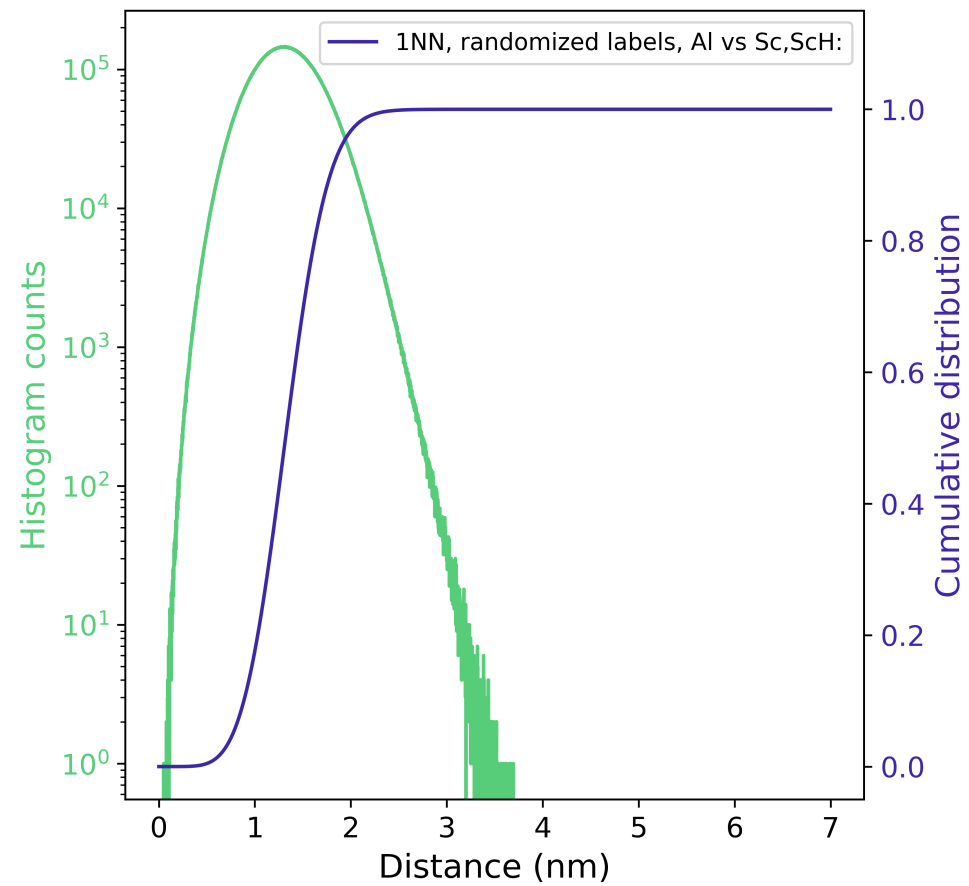

(a)


Figure 6: What is the kth nearest neighbor distribution function for particular ion types about specific ion types?



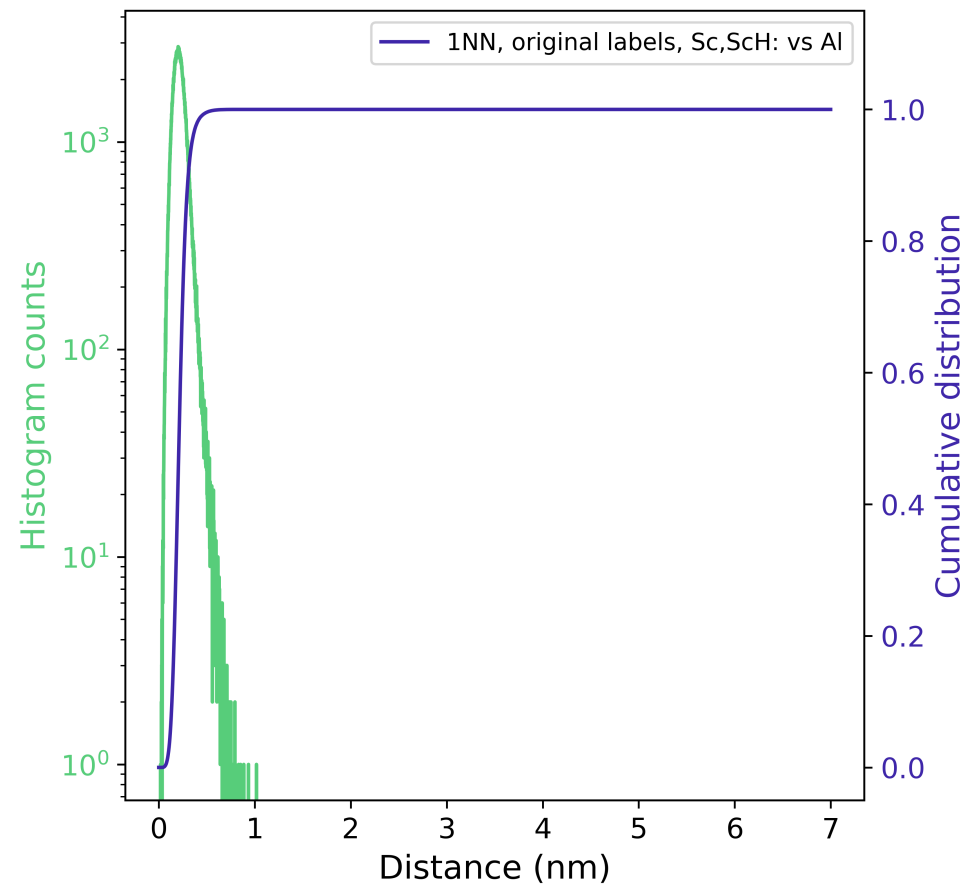


(a)

Figure 7: What is the kth nearest neighbor distribution function for particular ion types about specific ion types?

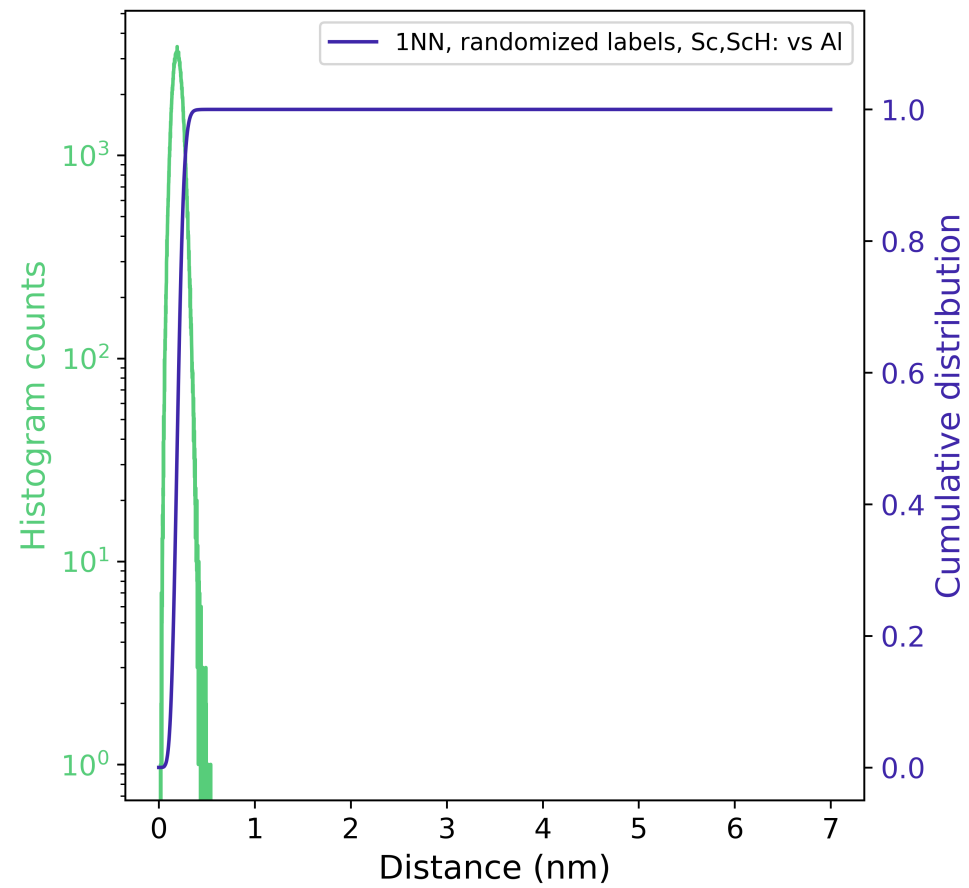


(a)

Figure 8: What is the kth nearest neighbor distribution function for particular ion types about specific ion types?

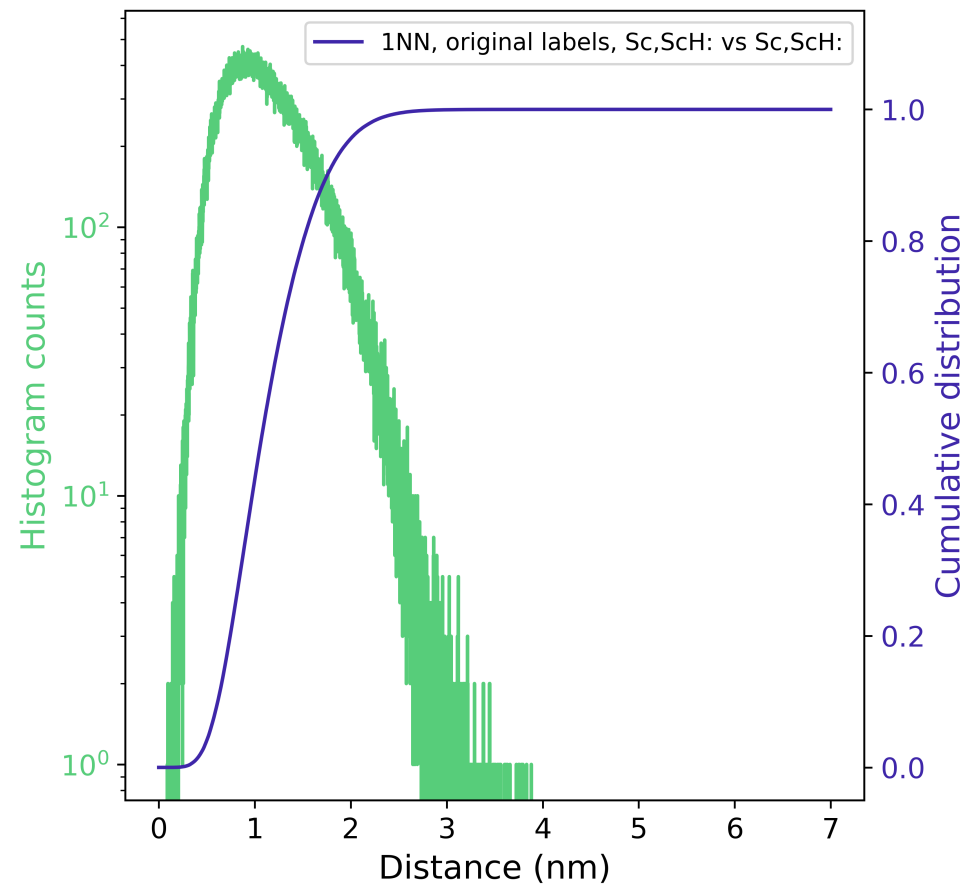


(a)

Figure 9: What is the kth nearest neighbor distribution function for particular ion types about specific ion types?

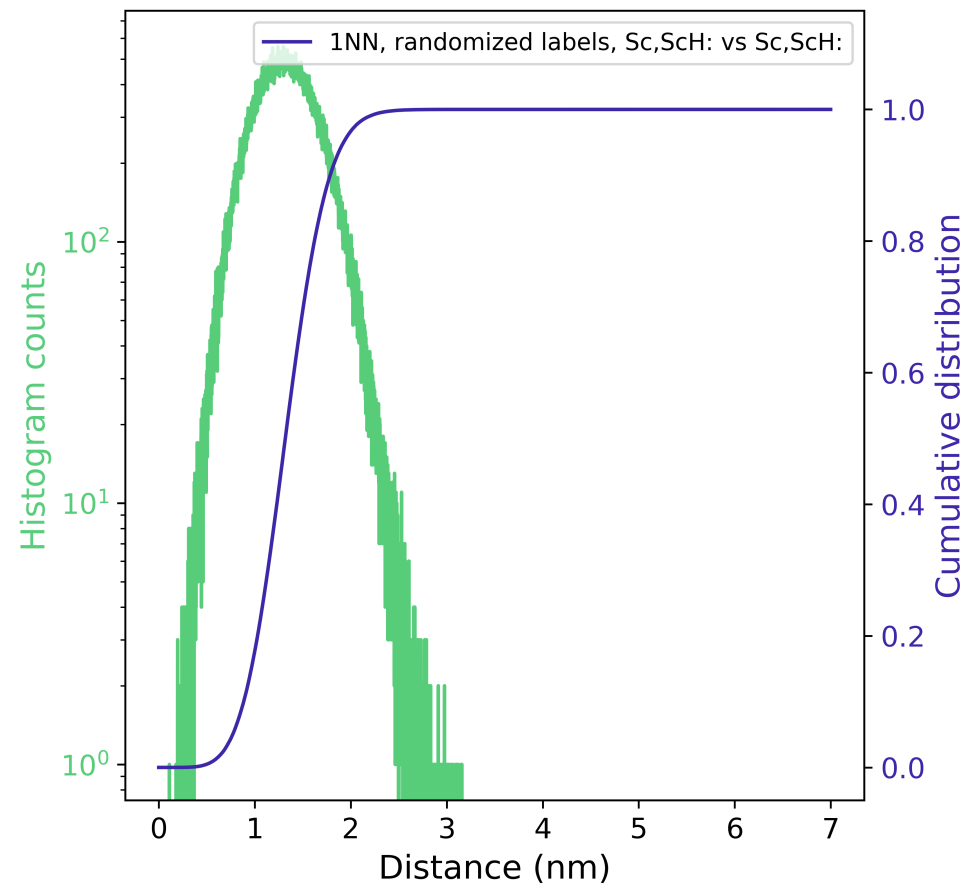


(a)

Figure 10: What is the kth nearest neighbor distribution function for particular ion types about specific ion types?

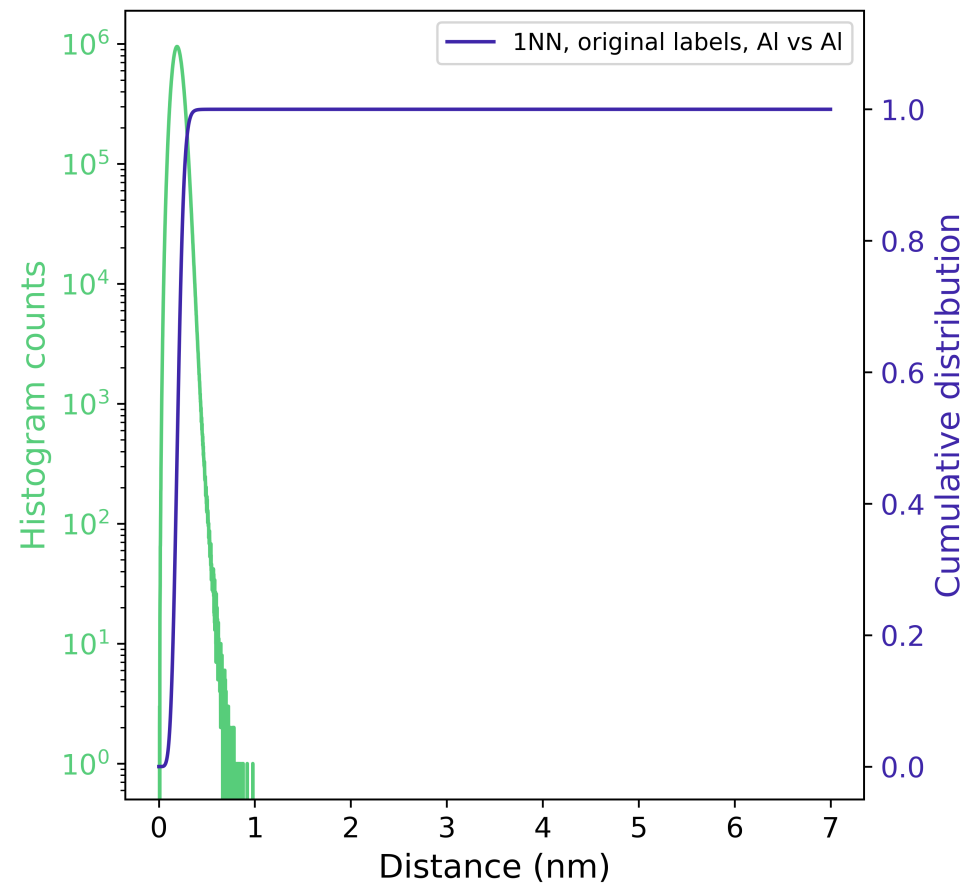

(a)

Figure 11: What is the kth nearest neighbor distribution function for particular ion types about specific ion types?

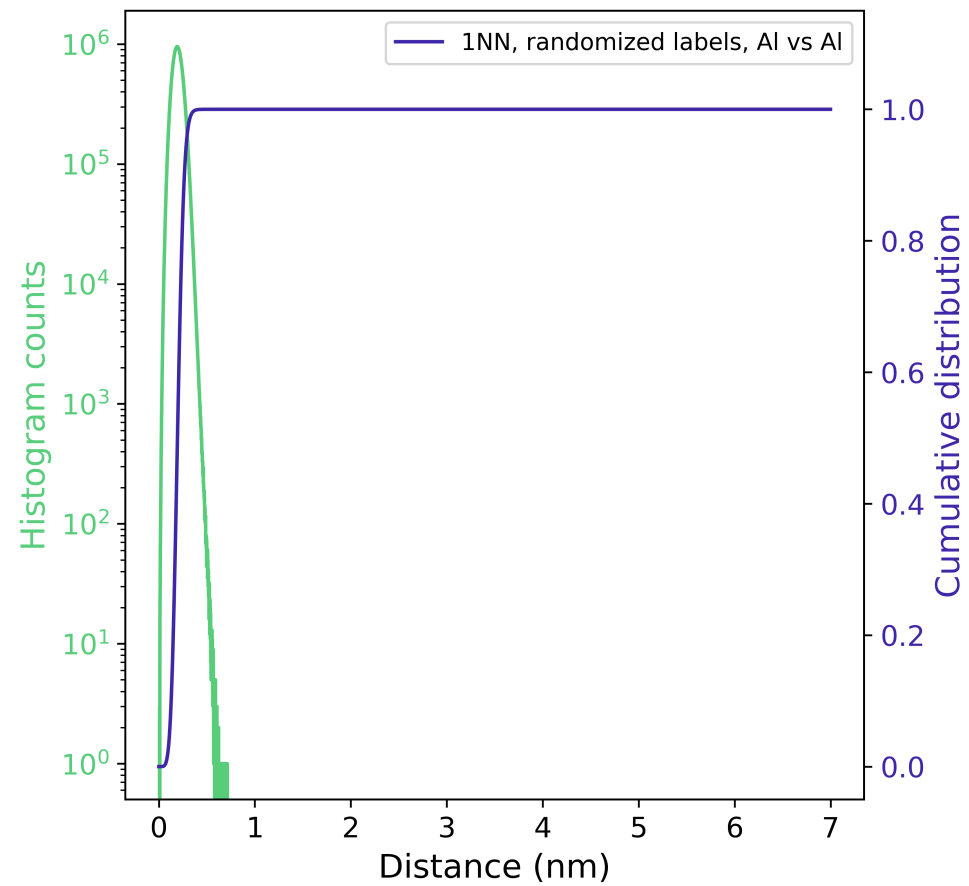


(a)

Figure 12: What is the kth nearest neighbor distribution function for particular ion types about specific ion types?

Table 6: Which XML settings were used for paraprobe-spatstat?

<b>Keyword</b>	<b>Value</b>	<b>Unit</b>	<b>Description</b>
AnalyzeKNN	1		
AnalyzeRDF	1		
AnalyzeRandomize	1		
AnalyzeSDM	0		
InputfileHullAndDistances	PARAPROBE.Surfacer.Results.SimID.3014001.h5		
InputfilePSE	PARAPROBE.PeriodicTableOfElements.xml		
InputfileReconstruction	PARAPROBE.Transcoder.Results.SimID.3014001.h5		
IontypeCombi0	Targets;Al;Neighbors;Al		
IontypeCombi1	Targets;Al;Neighbors;Sc,ScH:		
IontypeCombi2	Targets;Sc,ScH;;Neighbors;Al		
IontypeCombi3	Targets;Sc,ScH;;Neighbors;Sc,ScH:		
KOrderForKNN	1		
KOrderForSDM			
MaxSizeCachedResPerNode	17179869184		
PRNGType	MT19937		
PRNGWarmup	700000		
PRNGWorldSeed	18446744073697205938		
ROIRadiiKNNIncr	0.001	nm	
ROIRadiiKNNMax	7	nm	
ROIRadiiKNNMin	0	nm	
ROIRadiiRDFIncr	0.001	nm	
ROIRadiiRDFMax	7	nm	
ROIRadiiRDFMin	0	nm	
ROIRadiiSDMIncr	1	nm	
ROIRadiiSDMMax	0	nm	
ROIRadiiSDMMin	0	nm	
ROIVolumeInsideOnly	1		

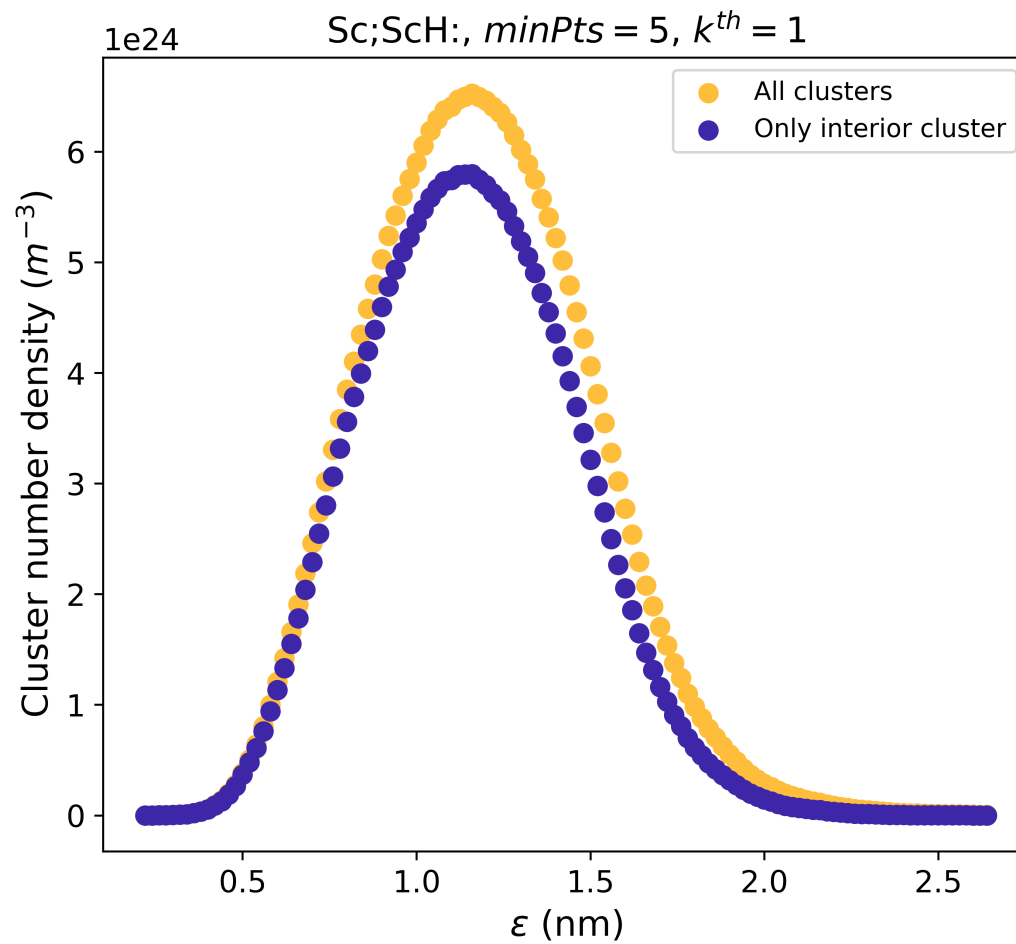

(a)

Figure 13: What are the number of clusters per unit volume, quantified using high-throughput studies of the DBScan/maximum separation clustering parameter?



Table 7: Which XML settings were used for paraprobe-dbscan?

Keyword	Value	Unit	Description
ClusteringMethod	2		
DBScanEpsilonIncr	0.02	nm	
DBScanEpsilonMax	7	nm	
DBScanEpsilonMin	0.2	nm	
DBScanMinPtsIncr	1		
DBScanMinPtsMax	1		
DBScanMinPtsMin	1		
DatasetEdgeThresholdDistance	1	nm	
IOStoreClusterIDs	0		
IOStoreClusters	1		
InputfileHullAndDistances	PARAPROBE.Surfacr.Results.SimID.3014001.h5		
InputfilePSE	PARAPROBE.PeriodicTableOfElements.xml		
InputfileReconstruction	PARAPROBE.Transcoder.Results.SimID.3014001.h5		
IontypeCombi0	Targets;Sc,ScH:		
MaxSepNumberOfIonsIncr	10		
MaxSepNumberOfIonsMax	5		
MaxSepNumberOfIonsMin	5		
MaxSizeCachedResPerNode	17179869184		

Table 8: How many MPI processes times OpenMP threads respectively were used, and what was the elapsed time for each tool run?

Toolname	#MPI	#OMP	$t_{elapsed}^{max}$ (s)
Transcoder	1	1	insignificant
Ranger	1	1	insignificant
Surfacr	8	40	36877.0
Tessellator	1	40	196.0
Spatstat	1	40	22120.0
DBScan	1	40	473.0